


ONE PIECE ORS FOR EVACTRON 45, 40, and 25 Models

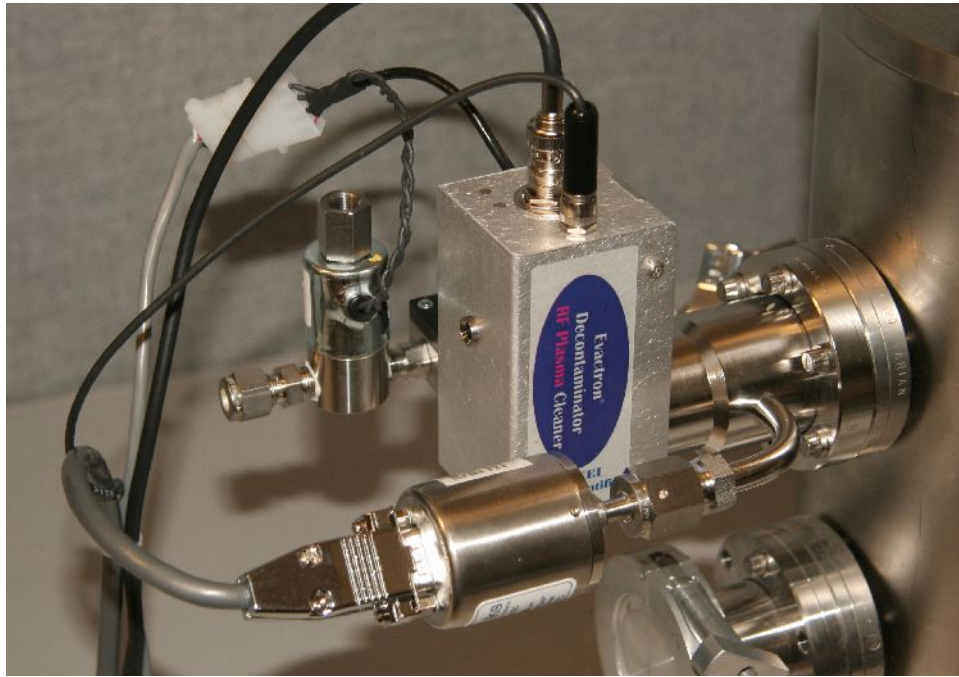

The One Piece ORS (Oxygen Radical Source) is an extra cost option for the Evactron 45 and 25 model series of Evactron Decontaminators. It has been designed for use where a cleaner look, fewer O-rings seals, and a compact space are requirements. Made of all stainless steel with VCR fittings for the gauges and valve it is suitable for operation to pressures down to 10^{-8} Torr in high vacuum equipment. All metal seals are used in the CF version except for viton valve seats in solenoid valves. VCR fittings are used for the Mini-Pirani gauge and to connect the valves. The external fitting is CF 2.75" or KF40.